

COZY 4.0

Far Infrared Throw/Wrap



**The Far Infrared Cozy™ 4.0 Wrap and Throw uses
Our *latest* most advanced Micron 4.0 Far infrared Technology.**

Using Your Cozy 4.0

1. Plug the connector into the Cozy 4.0 socket: 3 Steps

Step 1



Step 2



Step 3



2. Plug the controller mains cable with the 3 pin plug into the mains and then switch on the power at the mains plug. The display light will come on - The Cozy 4.0 is **not** heating yet.

Temperature Display



Power Button & Temperature control



Telephone: 01727 602022


Website: www.get-fitt.com

Email: admin@get-fitt.com

3. Turning on Your Cozy 4.0

Press and hold the selector button until a number shows on the display; then release the button.

Press & Hold the Power Button

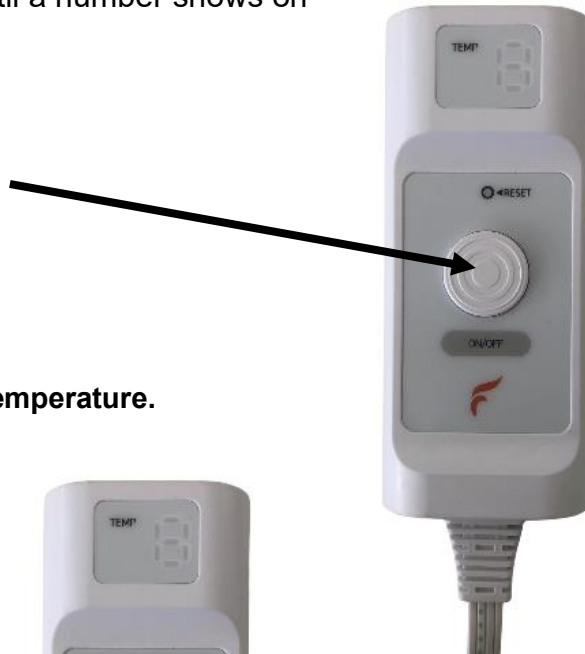

5. Press the selector to choose your temperature.

Press the Temperature Button to select your temperature



Once activated your Cozy 4.0 can take 2/3 minutes to reach the desired temperatures, this can take a little longer for the higher operating temperatures and may vary according to the ambient room temperature.

Temperature Settings

Setting	1	2	3	4	5	6	H
Temperature °C	35	40	45	50	55	60	65

Ending Your Session

When you wish to stop using the Cozy 4.0, press and hold the power button until the display shows a horizontal line. Then unplug when finished.

NOW unplug from the mains

Telephone: 01727 602022

Website: www.get-fitt.com

Email: admin@get-fitt.com